



Urea Agar Base (Christensen)(Autoclavable)

M112

Intended Use:

Urea Agar Base with the addition of Urea is recommended for the detection of urease production, particularly by members of the genus *Proteus* .

Composition**

Ingredients	Gms / Litre
Peptone	1.000
Dextrose (Glucose)	1.000
Sodium chloride	5.000
Disodium phosphate	1.200
Monopotassium phosphate	0.800
Phenol red	0.012
Agar	15.000
Final pH (at 25°C)	6.8±0.2

**Formula adjusted, standardized to suit performance parameters

Directions

Suspend 24.01 grams in 950 ml distilled water. Heat to boiling to dissolve the medium completely. Sterilize by autoclaving at 10 lbs pressure (115°C) for 20 minutes. Cool to 45-50°C and aseptically add 50 ml of sterile 40% Urea Solution (FD048) and mix well. Dispense into sterile tubes and allow to set in the slanting position. Do not overheat or reheat the medium as urea decomposes very easily.

Principle And Interpretation

Urea Agar is used to detect urease production. Urea Agar described by Christensen (3,7) detected urease activity by all rapidly urease-positive *Proteus* organisms and also by other members of *Enterobacteriaceae* (3) that exhibited a delayed urease reaction (8). This was accomplished by :

- adding glucose to the medium.
- decreasing the peptone concentration and
- decreasing the buffering system, as a less buffered medium detects even smaller amount of alkali (4).

Peptone is the source of essential nutrients. Dextrose is the energy source. Sodium chloride maintains the osmotic equilibrium of the medium whereas phosphates serve to buffer the medium. Urea is hydrolyzed to liberate ammonia. Phenol red indicator detects the alkalinity generated by visible colour change from orange to pink.

Prolonged incubation may cause alkaline reaction in the medium. A medium without urea serves as negative control to rule out false positive results. Also, all urea test media rely on the alkalinity formation and so they are not specific for determining the absolute rate of urease activity (8). The utilization of proteins may raise the pH to alkalinity due to protein hydrolysis and excess of amino acids liberation results in false positive reaction.

Type of specimen

Pure isolate from samples.

Specimen Collection and Handling

For clinical samples follow appropriate techniques for handling specimens as per established guidelines(5,6).

For food and dairy samples, follow appropriate techniques for sample collection and processing as per guidelines(1,9,10).

For water samples, follow appropriate techniques for sample collection, processing as per guidelines and local standards(2).

After use, contaminated materials must be sterilized by autoclaving before discarding.

Warning and Precautions

In Vitro diagnostic use only. Read the label before opening the container. Wear protective gloves/protective clothing/eye protection/face protection. Follow good microbiological lab practices while handling specimens and culture. Standard precautions as per established guidelines should be followed while handling clinical specimens. Safety guidelines may be referred in individual safety data sheets.

Limitations

1. Prolonged incubation may cause alkaline reaction in the medium.
2. Also, all urea test media rely on the alkalinity formation and so they are not specific for determining the absolute rate of urease activity (8).
3. The utilization of proteins may raise the pH to alkalinity due to protein hydrolysis and excess of amino acids liberation results in false positive reaction.

Performance and Evaluation

Performance of the medium is expected when used as per the direction on the label within the expiry period when stored at recommended temperature.

Quality Control

Appearance

Light yellow to light pink homogeneous free flowing powder

Gelling

Firm, comparable with 1.5% Agar gel

Colour and Clarity of prepared medium

Yellowish orange coloured clear to slightly opalescent gel forms in tubes as slants

Reaction

Reaction of 2.4% w/v aqueous solution at 25°C. pH : 6.8±0.2

pH

6.60-7.00

Cultural Response

Cultural characteristics observed on addition of sterile 40% Urea Solution (FD048) after an incubation at 35-37°C for 18-24 hours.

Cultural Response

Organism	Inoculum (CFU)	Urease
Cultural Response		
<i>Escherichia coli</i> ATCC 25922 (00013*)	50-100	negative reaction, no change
# <i>Klebsiella aerogenes</i> ATCC 13048 (00175*)	50-100	negative reaction, no change
<i>Klebsiella pneumoniae</i> ATCC 13883 (00097*)	50-100	positive reaction, cerise colour
<i>Proteus mirabilis</i> ATCC 25933	50-100	positive reaction, cerise colour
<i>Proteus vulgaris</i> ATCC 13315	50-100	positive reaction, cerise colour
<i>Salmonella Typhimurium</i> ATCC 14028 (00031*)	50-100	negative reaction, no change

Key : *Corresponding WDCM numbers.

Formerly known as *Enterobacter aerogenes*

Storage and Shelf Life

Store below 30°C in a tightly closed container and the prepared medium at 2 - 8°C. Use before expiry date on the label. On opening, product should be properly stored dry, after tightly capping the bottle in order to prevent lump formation due to the hygroscopic nature of the product. Improper storage of the product may lead to lump formation. Store in dry ventilated area protected from extremes of temperature and sources of ignition Seal the container tightly after use. Use before expiry date on the label. Product performance is best if used within stated expiry period.

Disposal

User must ensure safe disposal by autoclaving and/or incineration of used or unusable preparations of this product. Follow established laboratory procedures in disposing of infectious materials and material that comes into contact with clinical sample must be decontaminated and disposed of in accordance with current laboratory techniques (5,6).

Reference

1. American Public Health Association, Standard Methods for the Examination of Dairy Products, 1978, 14th Ed., Washington D.C.
2. Baird R.B., Eaton A.D., and Rice E.W., (Eds.), 2015, Standard Methods for the Examination of Water and Wastewater, 23rd ed., APHA, Washington, D.C.
3. Christensen W. B., 1946, J. Bacteriol., 52:461.
4. Farmer J. J. III, McWhorter A. C., Huntley G. A., Catignani J., J. Clin. Microbiol. 1975: 1 (1): 106-107.
5. Isenberg, H.D. Clinical Microbiology Procedures Handbook. 2nd Edition.
6. Jorgensen, J.H., Pfaller, M.A., Carroll, K.C., Funke, G., Landry, M.L., Richter, S.S and Warnock, D.W. (2015) Manual of Clinical Microbiology, 11th Edition. Vol. 1.
7. MacFaddin J. F., 1985, Media for Isolation-Cultivation-Identification-Maintenance of Medical Bacteria, Vol. 1, Williams and Wilkins, Baltimore, Md.
8. MacFaddin J. F., 2000, Biochemical Tests for Identification of Medical Bacteria, 3rd Ed., Williams and Wilkins, Baltimore, Md.
9. Salfinger Y., and Tortorello M.L. Fifth (Ed.), 2015, Compendium of Methods for the Microbiological Examination of Foods, 5th Ed., American Public Health Association, Washington, D.C.
10. Wehr H. M. and Frank J. H., 2004, Standard Methods for the Microbiological Examination of Dairy Products, 17th Ed., APHA Inc., Washington, D.C.


Revision : 03 / 2018



In vitro diagnostic medical device



CE Marking

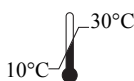

Storage temperature



Do not use if package is damaged



HiMedia Laboratories Pvt. Limited,
23 Vadhani Industrial Estate,
LBS Marg, Mumbai-86, MS, India



CE Partner 4U, Esdoornlaan 13, 3951
DB Maarn The Netherlands,
www.cepartner4u.eu

Disclaimer :

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory, diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.