



Semisolid Rappaport Vassiliadis Medium, Modified

M1282

Intended Use

Recommended for the detection of motile *Salmonella* species from food, faeces and environmental specimens.

Composition**

Ingredients	Gms / Litre
Tryptose	4.590
Tryptone	4.590
Sodium chloride	7.340
Potassium dihydrogen phosphate	1.470
Magnesium chloride, anhydrous	10.930
Malachite green oxalate	0.037
Agar	2.700
Final pH (at 25°C)	5.2±0.2

**Formula adjusted, standardized to suit performance parameters

Directions

Suspend 31.66 grams in 1000 ml purified / distilled water. Heat to boiling with frequent agitation to dissolve the medium completely. DO NOT AUTOCLAVE. Cool to 45-50°C and aseptically add the rehydrated contents of 1 vial of Novobiocin Selective Supplement (FD290) or 4 vials of Novobiocin supplement (FD096). Mix well and pour into sterile Petri plates. Air dry the plated medium at room temperature for at least one hour.

Principle And Interpretation

Semisolid Rappaport Vassiliadis Medium, Modified is prepared as per the formulation described by DeSmedt et al (2), for the detection of motile *Salmonella* species from food and environmental specimens. This medium detects more *Salmonella* positive samples than the routinely used enrichment procedures (3,4,5).

Tryptose, tryptone provide the nitrogenous and carbonaceous substances and other essential growth nutrients. Phosphate gives good buffering capacity to the medium. Malachite green oxalate and Novobiocin inhibits many gram positive bacteria. The working of medium is based on the ability of *Salmonella* species to migrate in the selective medium by competing with the other motile organisms, thus producing opaque halos of growth. This medium can be used in combination with direct culture and Selenite F Broth (M052A) enrichment for isolation of *Salmonella* species from faeces and sub-culturing on XLD Agar (M031) or Mannitol Lysine Agar (M1071) results in higher recovery rates (1). This medium is not suitable for the detection of non-motile strains of *Salmonella* (6).

Inoculate 3 drops (0.1 ml) of pre-enrichment culture (16-20 hours old) in separate spots on the air-dried medium surface. Incubate the plates in an upright position at 42°C for upto 24 hours. The motile bacteria will show a halo or zone of growth originating from inoculation spot. *Salmonella* species show straw coloured colonies. Sub-cultures can be carried out from the outside edge of the halo to confirm purity and for further biochemical and serological tests.

Type of specimen

Clinical samples - faeces, Food samples: meat and meat products.

Specimen Collection and Handling:

For clinical samples follow appropriate techniques for handling specimens as per established guidelines (7,8).

For food samples, follow appropriate techniques for sample collection and processing as per guidelines (9).

After use, contaminated materials must be sterilized by autoclaving before discarding.

Warning and Precautions

Read the label before opening the container. Wear protective gloves/protective clothing/eye protection/ face protection. Follow good microbiological lab practices while handling specimens and culture. Standard precautions as per established guidelines should be followed while handling specimens. Safety guidelines may be referred in individual safety data sheets.

Limitations :

1. Due to variable nutritional requirements, some strains show poor growth on this medium.
2. The plates should be incubated in upright and steady position as medium contains less agar.
3. Air dried plate should be used, or it may give erroneous results.

Performance and Evaluation

Performance of the medium is expected when used as per the direction on the label within the expiry period when stored at recommended temperature.

Quality Control

Appearance

Light yellow to light blue homogeneous free flowing powder

Gelling

Semisolid, comparable with 0.27% Agar gel.

Colour and Clarity of prepared medium

Greenish blue coloured clear to slightly opalescent semisolid gel forms in Petri plates which may have a slight precipitate.

Reaction

Reaction of 3.16% w/v aqueous solution at 25°C. pH : 5.2±0.2

pH

5.00-5.40

Cultural Response

M1282: Cultural characteristics observed after an incubation at 42-43°C for 18-24 hours with added Novobiocin Selective Supplement (FD290) or Novobiocin supplement (FD096).

Organism	Inoculum (CFU)	Growth at 42±1°C	Motility
<i>Citrobacter freundii</i> ATCC 8090 (00006*)	50-100	none-poor	negative reaction, no colour change
<i>Escherichia coli</i> ATCC 25922 (00013)*	50-100	none-poor	negative reaction, no colour change
<i>Salmonella Typhi</i> ATCC 6539	50-100	good-luxuriant	positive reaction, colourless to light pink zone
<i>Salmonella Typhimurium</i> ATCC 14028 (00031)*	50-100	good-luxuriant	positive reaction, colourless to light pink zone
<i>Salmonella Enteritidis</i> ATCC 13076 (00030*)	50-100	good-luxuriant	positive reaction, colourless to light pink zone

Key : *Corresponding WDCM numbers.

Storage and Shelf Life

Store between 10-30°C in a tightly closed container and the prepared medium at 2-8°C. Use before expiry date on the label. On opening, product should be properly stored dry, after tightly capping the bottle in order to prevent lump formation due to the hygroscopic nature of the product. Improper storage of the product may lead to lump formation. Store in dry ventilated area protected from extremes of temperature and sources of ignition. Seal the container tightly after use. Use before expiry date on the label. Product performance is best if used within stated expiry period.

Disposal

User must ensure safe disposal by autoclaving and/or incineration of used or unusable preparations of this product. Follow established laboratory procedures in disposing of infectious materials and material that comes into contact with sample must be decontaminated and disposed of in accordance with current laboratory techniques (7, 8).

Reference

1. Aspinall S.T., Hindle M.A. and Hutchinson D.N., 1992, Europ. J. Clin. Microbiol. Inf. Dis., 11:936.
2. De Smedt J.M., Balderdijk R., Rappold H. and Lautenschlaeger D., 1986, J. Food Prot., 49:510.
3. De Smedt J.M., Balderdijk R., 1987, J. Food Prof., 50:658.
4. De Zutter L. et al, 1991, Int. J. Food Microbiol., 13:11.
5. De Smedt J.M. et al, 1991, Int. J. Food Microbiol., 13:301.
6. Holbrook R. et al, 1989, Lett. Appl. Microbiol., 8:139.
7. Isenberg, H.D. Clinical Microbiology Procedures Handbook. 2nd Edition.
8. Jorgensen, J.H., Pfaller, M.A., Carroll, K.C., Funke, G., Landry, M.L., Richter, S.S and Warnock, D.W. (2015) Manual of Clinical Microbiology, 11th Edition. Vol. 1.
9. Salfinger Y., and Tortorello M.L. Fifth (Ed.), 2015, Compendium of Methods for the Microbiological Examination of Foods, 5th Ed., American Public Health Association, Washington, D.C.


Revision : 04 / 2018



In vitro diagnostic medical device



CE Marking

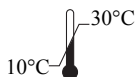

Storage temperature



Do not use if package is damaged



HiMedia Laboratories Pvt. Limited,
23 Vadhani Industrial Estate,
LBS Marg, Mumbai-86, MS, India



CE Partner 4U, Esdoornlaan 13, 3951
DB Maarn The Netherlands,
www.cepartner4u.eu

Disclaimer :

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory, diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.