



Anaerobic Agar

M228

Intended Use:

Anaerobic Agar is recommended for the cultivation of anaerobic bacteria, especially *Clostridium* species and other anaerobic organisms from clinical and non-clinical samples.

Composition**

Ingredients	Gms / Litre
Tryptone	20.000
Dextrose (Glucose)	10.000
Sodium chloride	5.000
Sodium thioglycollate	2.000
Sodium formaldehyde Sulfoxylate	1.000
Methylene blue	0.002
Agar	20.000
Final pH (at 25°C)	7.2±0.2

**Formula adjusted, standardized to suit performance parameters

Directions

Suspend 58.0 grams in 1000 ml distilled water. Heat to boiling to dissolve the medium completely. Sterilize by autoclaving at 15 lbs pressure (121°C) for 15 minutes. Cool to 45-50°C. Mix well and pour into sterile Petri plates.

Principle And Interpretation

Anaerobic Agar was originally designed for surface cultivation of members of the genus *Clostridium* and other anaerobic organisms on plates (1). This medium is suitable for isolation of facultative and obligate anaerobes and for the study of colonial morphology as colonies can be readily seen on the light coloured agar and are easily accessible (2,3). Anaerobic bacteria vary in their sensitivity to oxygen and nutritional requirements (3). Anaerobic bacteria lack cytochromes and thus are unable to use oxygen as a terminal electron acceptor (4).

This medium contains sodium thioglycollate and sodium formaldehyde sulfoxylate that provide adequate anaerobiosis which is indicated by methylene blue present in the medium which yields blue colour to medium in presence of oxygen. Casein enzymic hydrolysate and dextrose provide essential nutrients while sodium chloride maintains osmotic equilibrium.

Dispense 50-60 ml medium per 95 x 20 mm plate. For best results, use porous tops for the plates during solidification to get the dry surface. Inoculation can be done by streaking or smearing. Cover the inoculated plate with sterile Brewer Anaerobic Petri dish cover. Incubate aerobically, as desired. When standard plates are used, dispense 0.1 to 1.0 ml of inoculum into plates and mix with 20 - 25 ml of sterile medium. After solidification, incubate anaerobically as required by particular organism under study. Methylene blue is inhibitory to some anaerobic microorganisms.

Type of specimen

Clinical- stool, blood, abscess

Specimen Collection and Handling

For clinical samples follow appropriate techniques for handling specimens as per established guidelines (3,5).

After use, contaminated materials must be sterilized by autoclaving before discarding.

Warning and Precautions

In Vitro diagnostic use only. Read the label before opening the container. Wear protective gloves/protective clothing/ eye protection/face protection. Follow good microbiological lab practices while handling specimens and culture. Standard precautions as per established guidelines should be followed while handling clinical specimens. Safety guidelines may be referred in individual safety data sheets.

Please refer disclaimer Overleaf.

Limitations

- 1.Ensure that the clinical samples are properly transported under anaerobic conditions.
- 2.Proper anaerobic conditions must be maintained for optimal recovery of organisms
3. Methylene blue is toxic to certain anaerobes

Performance and Evaluation

Performance of the medium is expected when used as per the direction on the label within the expiry period when stored at recommended temperature.

Quality Control

Appearance

Cream to yellow homogeneous free flowing powder

Gelling

Firm, comparable with 2.0% agar gel.

Colour and Clarity of prepared medium

Light amber coloured, clear to slightly opalescent gel forms in Petri plates that becomes greenish due to aeration on standing

Reaction

Reaction of 5.8% w/v aqueous solution at 25°C. pH : 7.2±0.2

pH

7.00-7.40

Cultural Response

Cultural characteristics observed under anaerobic condition after an incubation at 35-37°C for 48-72 hours.

Organism	Inoculum (CFU)	Growth	Recovery
Cultural Response			
<i>Clostridium perfringens</i> ATCC 12924	50-100	good-luxuriant	≥50%
<i>Clostridium sporogenes</i> ATCC 11437	50-100	good-luxuriant	≥50%
<i>Clostridium butyricum</i> ATCC 13732	50-100	good-luxuriant	≥50%

Storage and Shelf Life

Store below 30°C in a tightly closed container and the prepared medium at 2 - 8°C. Use before expiry date on the label. On opening, product should be properly stored dry, after tightly capping the bottle in order to prevent lump formation due to the hygroscopic nature of the product. Improper storage of the product may lead to lump formation. Store in dry ventilated area protected from extremes of temperature and sources of ignition Seal the container tightly after use. Use before expiry date on the label. Product performance is best if used within stated expiry period.

Disposal

User must ensure safe disposal by autoclaving and/or incineration of used or unusable preparations of this product. Follow established laboratory procedures in disposing of infectious materials and material that comes into contact with clinical sample must be decontaminated and disposed of in accordance with current laboratory techniques (3,5).

Reference


1. Brewer J. H., 1942, Science, 95:587.
2. Vera J., 1942, J. Bacteriol., 44:497.
3. Isenberg (Ed.), 1992, Clinical Microbiology Procedures Handbook, American Society for Microbiology, Washington, D.C.
4. Baron E. J., Peterson and Finegold S. M., Bailey & Scotts Diagnostic Microbiology, 9th Ed., 1994, Mosby-Year Book, Inc., St. Louis, Mo.
5. Jorgensen, J.H., Pfaller, M.A., Carroll, K.C., Funke, G., Landry, M.L., Richter, S.S and Warnock., D.W. (2015) Manual of Clinical Microbiology, 11th Edition. Vol. 1.



In vitro diagnostic medical device



CE Marking

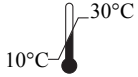

Storage temperature



Do not use if package is damaged



HiMedia Laboratories Pvt. Limited,
B /4-6 , MIDC, Dindori, Nashik MH

www.himedialabs.com



CE Partner 4U ,Esdoornlaan 13, 3951
DB Maarn The Netherlands,
www.cepartner4u.eu

Disclaimer :

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory, diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.